





IonOptix FHD Cell Stimulation/Perfusion Chamber

Bath Width:
6 mm

Bath Length:
19 mm

Bath Depth:
3.5 mm

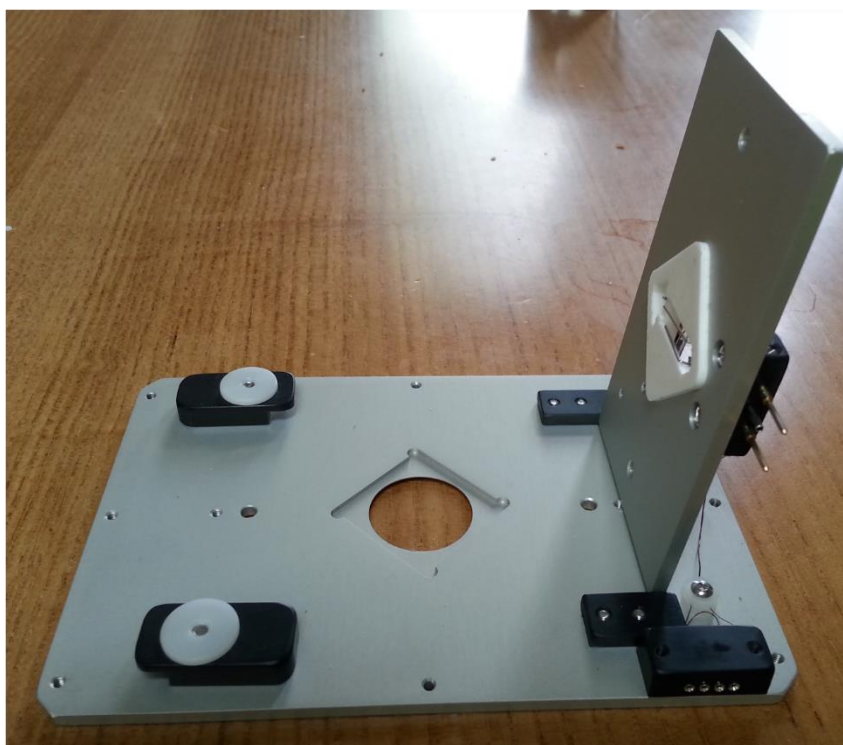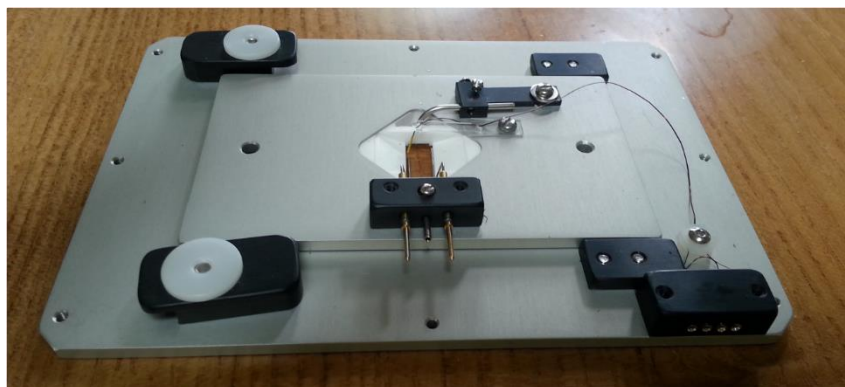




enabling discovery

Platinum lead wires

Perfusion outlet
w/ adj. height to
change bath volume


10kohm thermistor

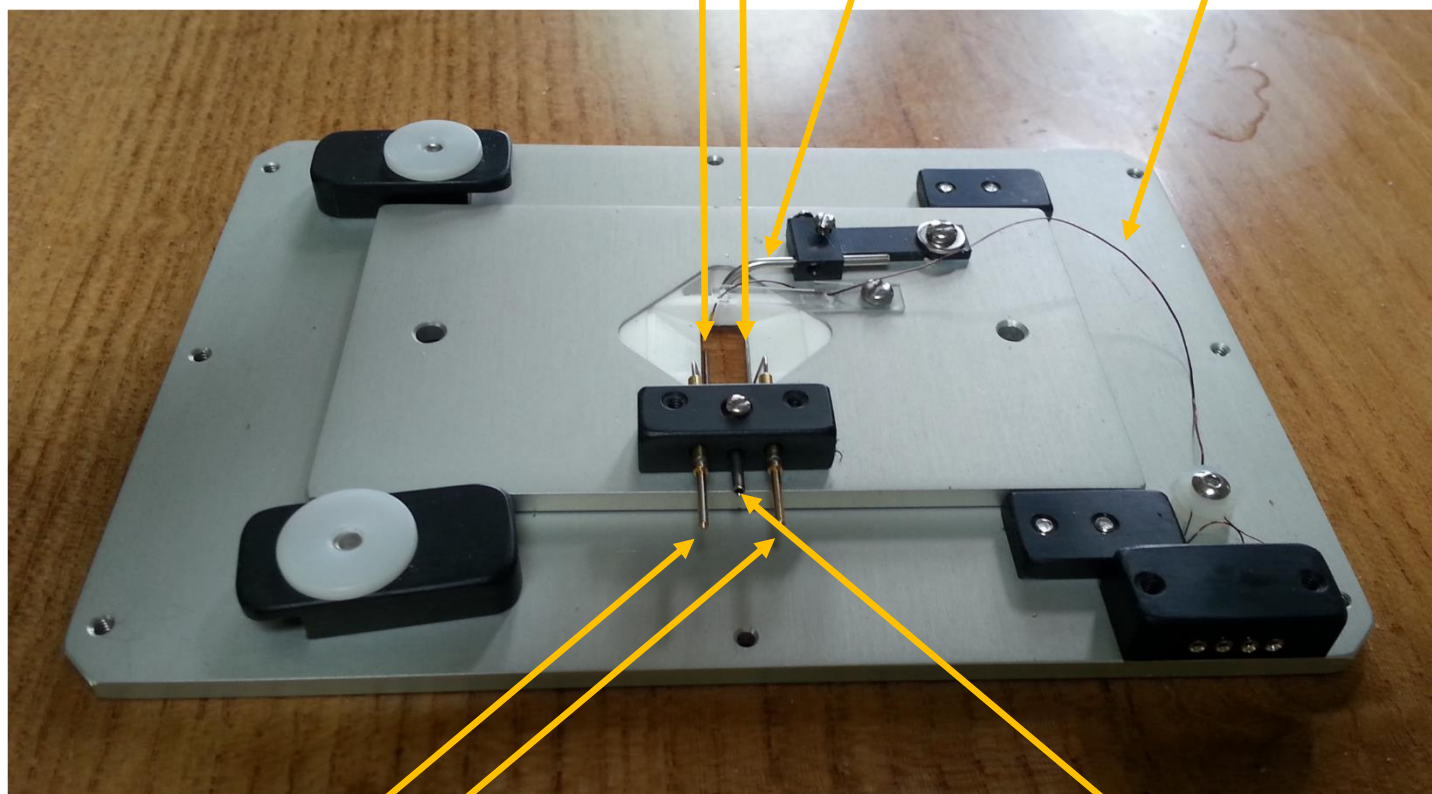

1 mm pin connectors

Perfusion inlet

Gasket seals bath and provides constrained geometry for laminar flow

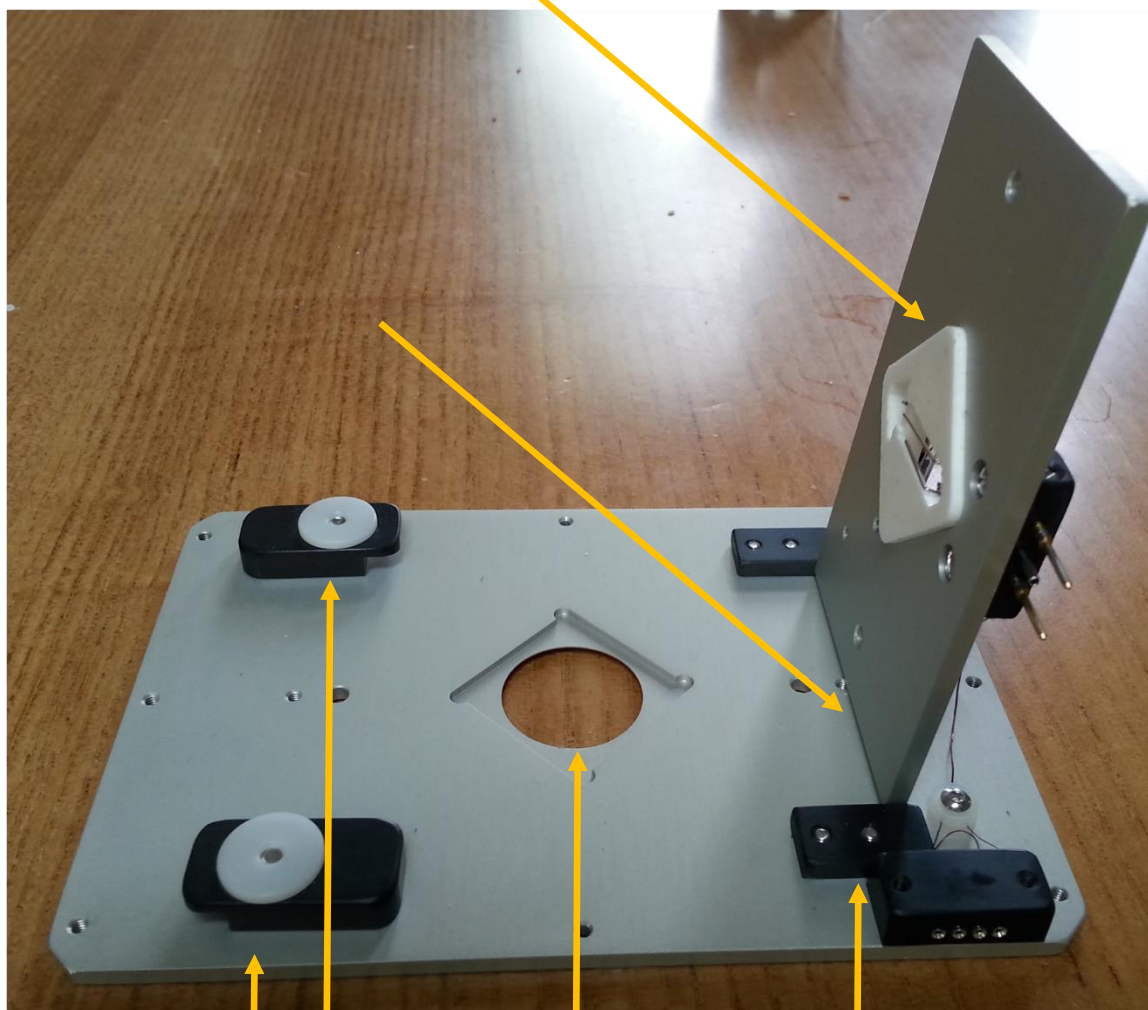

Lid Locks twist to close lid

Well for 25 mm cover glass

Hinged lid for rapid replacement of
cover glass and cleaning